

USING YOUR NEW LOOP EARPLUGS

Quiet | Quiet Plus | Experience | Experience Plus | Engage | Engage Plus | Engage Kids | Dream | Switch



Quiet

Quiet Plus

Engage

Engage Plus

Experience

Experience Plus

Switch

How to

INSERT YOUR EARPLUGS

- Insert the earplug with two fingers
- Twist to the back until it sits at the base of your ear comfortably
- The earplugs shouldn't stick out or be visible from the front




All earplugs


How to

FIND YOUR PERFECT SIZE


Quiet | Engage | Experience | Engage Kids | Switch

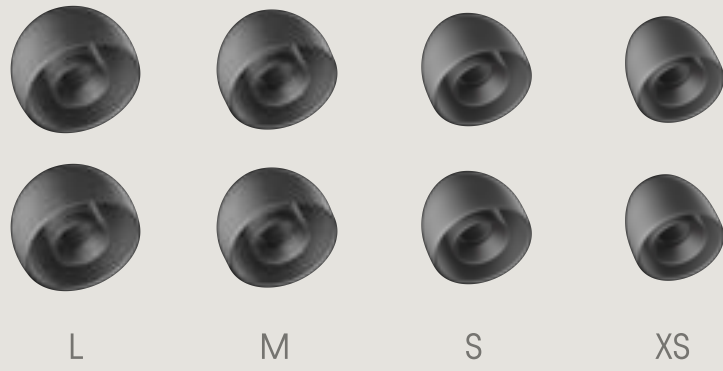


Engage Plus | Experience Plus

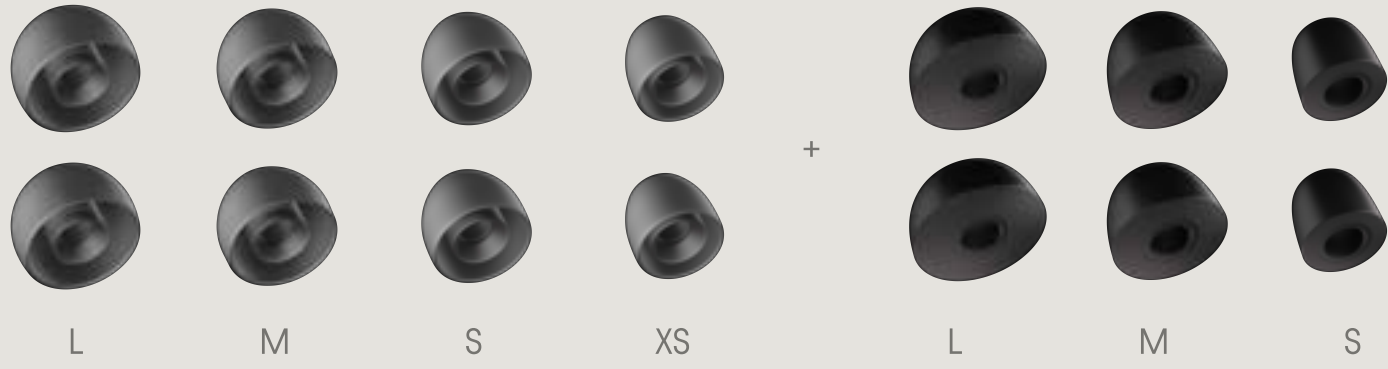

Quiet Plus

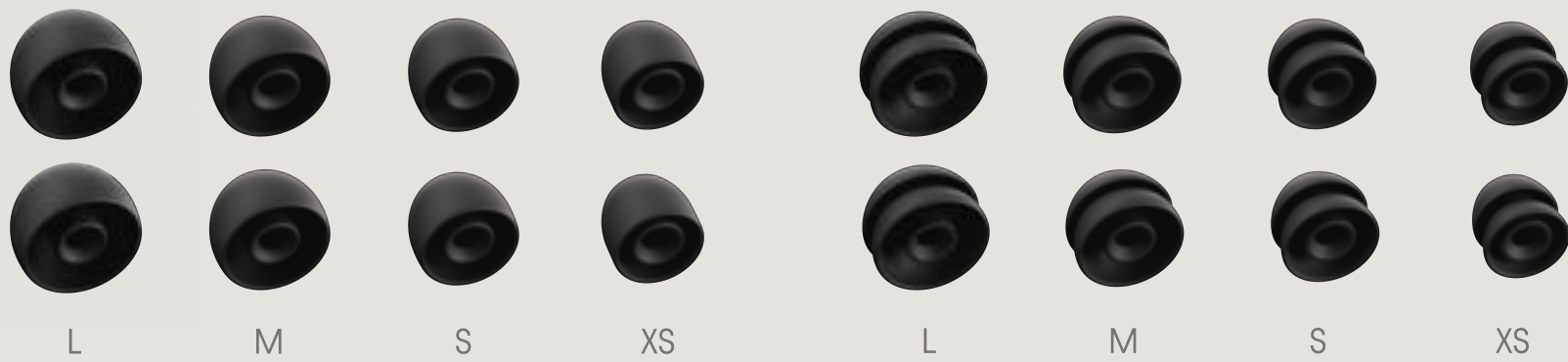


Dream

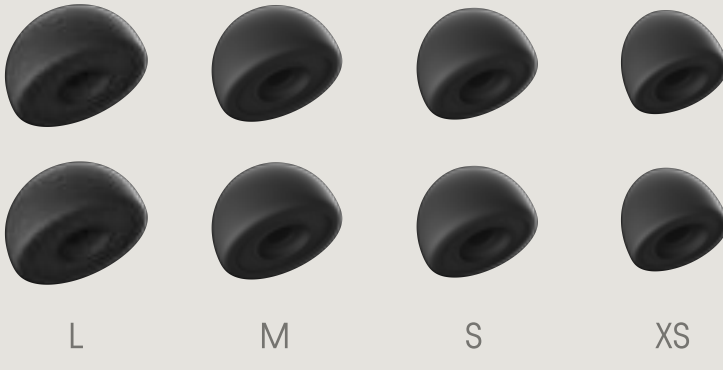

All earplugs

- > Try all ear tip sizes to find your perfect fit
- > If they don't feel secure, try a size up
- > If they stick out or cause discomfort, try a smaller size

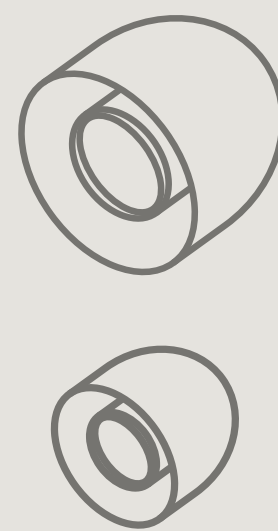

Left and right tip sizes can differ

Hearing yourself talk more loudly
is a sign of a good seal




Loop Switch

How to SWITCH MODES

Slide the dial up or down while Loop Switch is in your ear to select your preferred noise reduction mode.

Choose Quiet Mode for maximum noise reduction, Engage Mode for socializing and Experience Mode for live music and events.

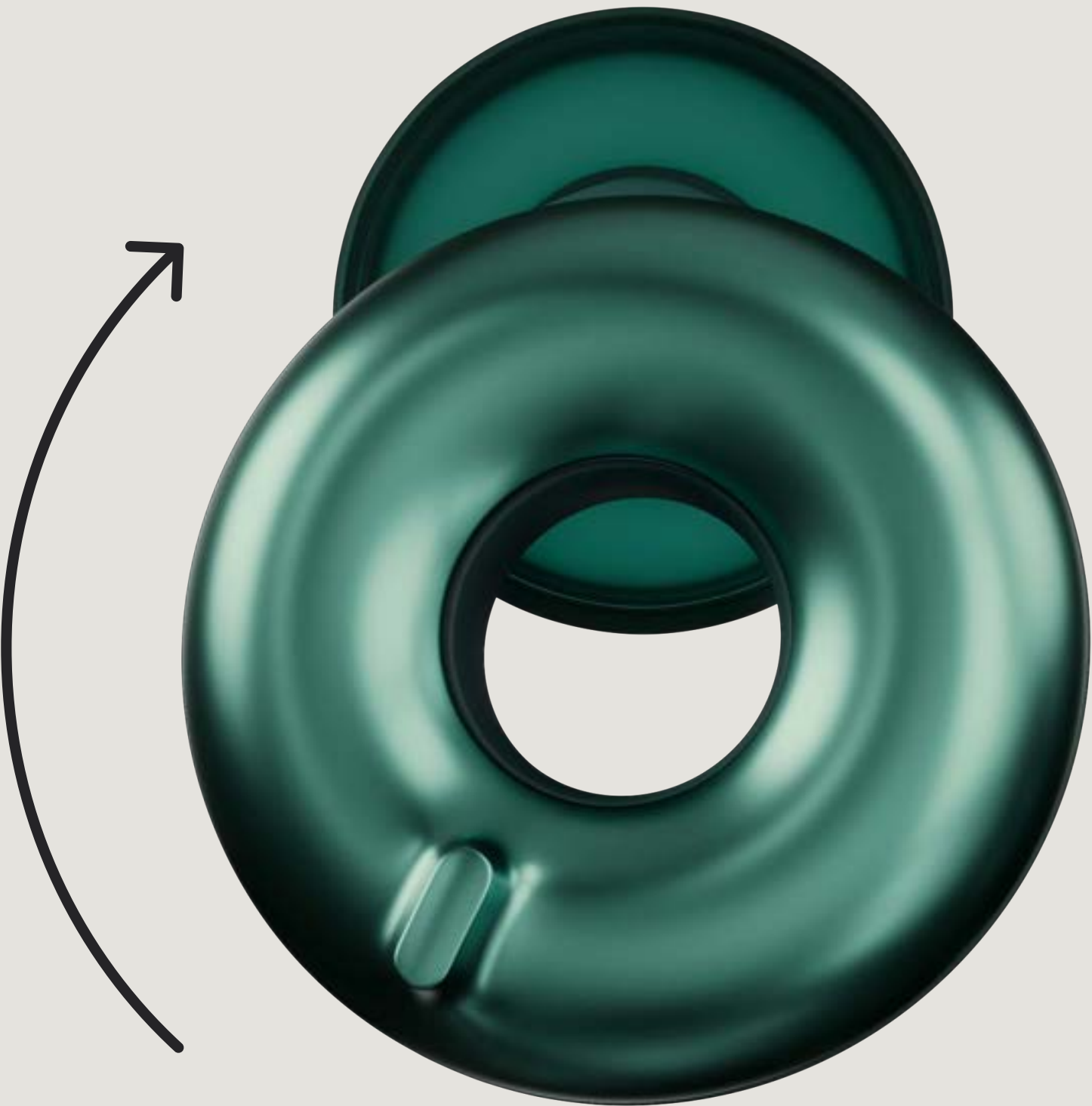

	● Quiet Mode	◐ Experience Mode	○ Engage Mode
NRR	20 dB	17 dB	24 dB
SNR	26 dB	23 dB	20 dB